




PERIPHERAL INTERVENTION

ONE-STEP™

Centesis Catheter

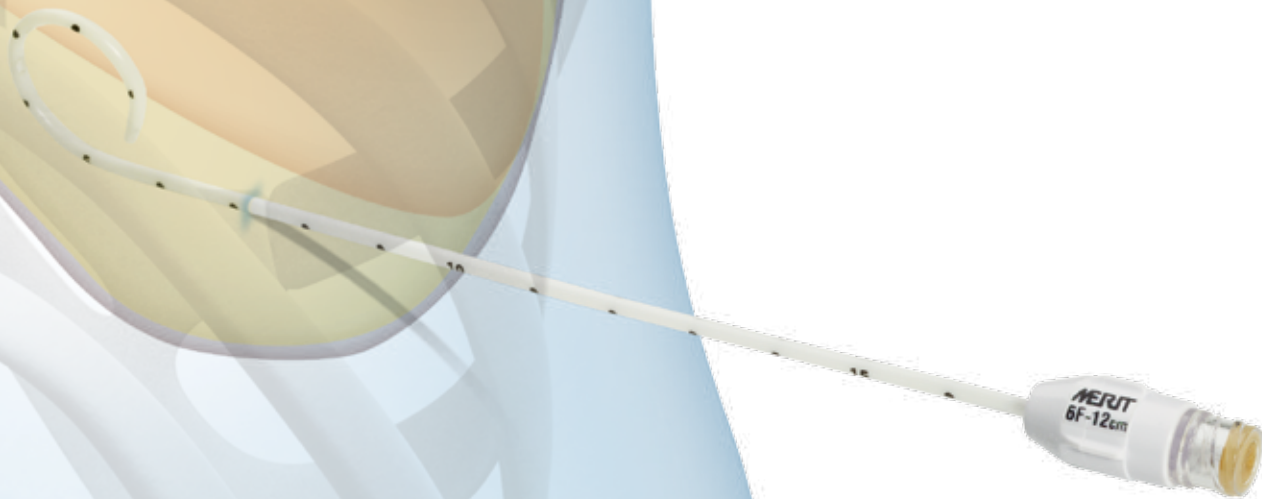

Up to **34% LESS** catheter
tip penetration force¹ and up
to **16% FASTER** flow rate^{1,2}



ONE-STEP™

Centesis Catheter

Safe, controlled insertion with decreased procedure times.

The valved One-Step centesis catheter is designed with an integrated self-sealing valve to minimize the risk of air entering the pleural space.



Pigtail shape
reduces drainage hole occlusion

Distal obturator
prevents catheter puncture and shearing and accidental needlestick

Large holes
increase flow rates to decrease procedure times

Depth markers and echo-enhanced needle
for fast and accurate positioning



Self-closing valve
prevents air infiltration and fluid leakage

Kink-resistant material
keeps catheter open to help maximize drainage

Sharp needle and smooth transitions
for ease of insertion and decreased patient discomfort



PERIPHERAL INTERVENTION

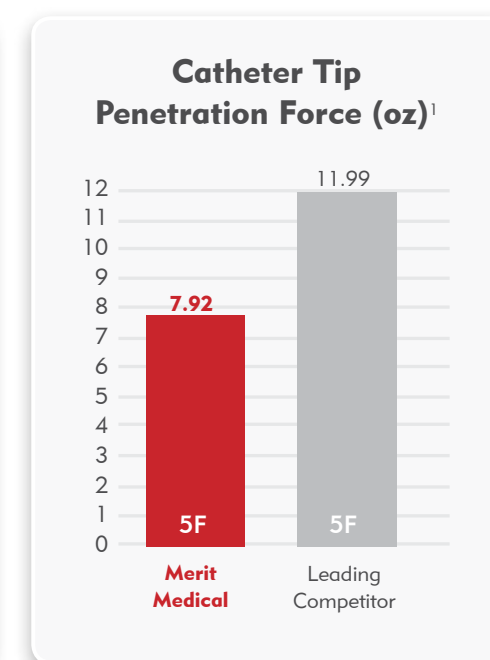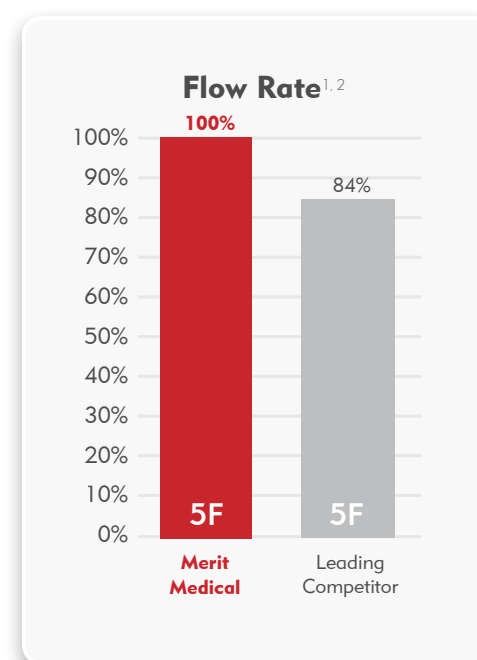

1. Testing as of 2015. Data on file.
2. Tested with 10 cm straight catheters.

ONE-STEP™

Centesis Catheter

TRAY MAY INCLUDE

| | |
|----|--|
| 1 | Valved One-Step Centesis Catheter |
| 1 | ChloroPrep® 3 mL tinted* |
| 1 | Drape tri-layer |
| 1 | Fenestrated drape |
| 1 | Lidocaine 20 mL |
| 1 | Lidocaine label |
| 1 | 18G x 1" shielded needle |
| 1 | 22G x 1.5" shielded needle |
| 1 | 25G x 1" shielded needle |
| 1 | Medallion® Syringe 10 mL |
| 1 | 10 mL Luer lock syringe |
| 1 | 60 mL Luer lock syringe |
| 1 | #11 Futura® Safety Scalpel |
| 1 | ShortStop® Temporary Sharps Holder |
| 10 | Gauze 4" x 4" (10 cm x 10 cm) |
| 1 | Non-vented female cap |
| 3 | Specimen vials (10 mL) with caps |
| 1 | 3-Way large bore stopcock |
| 1 | 48" (122 cm) tubing with spike for evacuated container |
| 1 | 96" (244 cm) wall suction tubing |
| 1 | Stepped adapter |
| 1 | Bandage dressing |
| 1 | CSR wrap |



*ChloroPrep® is a registered trademark of CareFusion 2200, Inc.

ORDERING INFORMATION

| Part Number | French Size | Length (cm) | Tip | Valve | Luer | Drainage Holes |
|----------------------|-------------|-------------|----------|-------|-------|----------------|
| CATHETER ONLY | | | | | | |
| OSC-5F-12PLT | 5 | 12 | Pigtail | Yes | Fixed | 4 |
| OSC-5F-12PL | 5 | 12 | Pigtail | - | Fixed | 4 |
| OSC-5F-7LT | 5 | 7 | Straight | Yes | Fixed | 4 |
| OSC-5F-10LT | 5 | 10 | Straight | Yes | Fixed | 4 |
| OSC-5F-7L | 5 | 7 | Straight | - | Fixed | 4 |
| OSC-5F-10L | 5 | 10 | Straight | - | Fixed | 4 |
| OSC-6F-12PLT | 6 | 12 | Pigtail | Yes | Fixed | 4 |
| OSC-6F-12PL | 6 | 12 | Pigtail | - | Fixed | 4 |
| OSC-6F-7LT | 6 | 7 | Straight | Yes | Fixed | 4 |
| OSC-6F-10LT | 6 | 10 | Straight | Yes | Fixed | 4 |
| TRAY | | | | | | |
| TAPS-5F-12PLT | 5 | 12 | Pigtail | Yes | Fixed | 4 |
| SPPT-5F-12PL | 5 | 12 | Pigtail | - | Fixed | 4 |
| TAPS-5F-10LT | 5 | 10 | Straight | Yes | Fixed | 4 |
| SPPT-5F-7L | 5 | 7 | Straight | - | Fixed | 4 |
| SPPT-5F-10L | 5 | 10 | Straight | - | Fixed | 4 |
| TAPS-6F-12PLT | 6 | 12 | Pigtail | Yes | Fixed | 4 |
| SPPT-6F-12PL | 6 | 12 | Pigtail | - | Fixed | 4 |
| TAPS-6F-10LT | 6 | 10 | Straight | Yes | Fixed | 4 |
| SPPT-6F-10L | 6 | 10 | Straight | - | Fixed | 4 |

Packaged 5 per box. Not all products listed. **4F and 8F also available.**
For a complete inventory please visit our website or talk to your Merit representative.



Understand. Innovate. Deliver.™

Merit Medical Systems, Inc.
1600 West Merit Parkway
South Jordan, Utah 84095
1.801.208.4300
1.800.35.MERIT

Merit Medical Europe, Middle East, & Africa (EMEA)
Amerikalaan 42, 6199 AE Maastricht-Airport
The Netherlands
+31 43 358 82 22

Merit Medical Ireland Ltd.
Parkmore Business Park West
Galway, Ireland
+353 (0) 91 703 733

Merit.com

Austria
0800 295 374

Belgium
0800 72 906 (Dutch)
0800 73 172 (Français)

Denmark
80 88 00 24

Finland
0800 770 586

France
0800 91 60 30

Germany
0800 182 0871

Ireland (Republic)
1800 553 163

Italy
800 897 005

Luxembourg
8002 25 22

Netherlands
0800 022 81 84

Norway
800 11629

Portugal
308 801 034

Russia
+7 495 221 89 02

Spain
+34 911238406

Sweden
020 792 445

Switzerland
(Deutsch)
+41 225180252
(Français)
+41 225948000
(Italiano)
+41 225180035
UK
0800 973 115